# ARDESTO

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## **IMPORTANT SAFEGUARDS**

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### **READ ALL INSTRUCTIONS BEFORE USING.**

#### **DANGER:**

To reduce the risk of electrocution:

Do not handle charger with wet hands.

Do not place charger in or drop charger into water or other liquid.

Check the charger cord for damage before the first use and during the life of the product.

Use this product only as indicated in these instructions or as recommended by your dental professional.

Only use tips and accessories that are recommended by manufacturer.

Do not drop or insert any foreign object into any opening on the product.

Keep charging unit and cord away from heated surfaces.

Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product.

Use of these can reduce performance and will shorten the life of the product. Use of these products may void warranty.

Fill reservoir with water or other dental professional recommended solution only.

Remove any oral jewelry prior to use of this product.

Do not use if you have an open wound on your tongue or in your mouth.

**WARNING:**

**To reduce the risk of burns, electrocution, fire, or injury to persons:**

Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger. The device must only be powered by a safety extra-low voltage as indicated on the label on the device.

Improper voltage will damage product and VOID warranty.

If product is opened/disassembled for ANY reason, warranty is VOID.

Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Please contact authorized service center for help. If the power cord is damaged, it must be replaced with a dedicated cord or kit provided by the manufacturer or service.

Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.

Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.

If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any

other oral hygiene aid.

This device can be used by children over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are under constant supervision or instructed in the safe use of the device and they understand the possible dangers. Cleaning and customer service should not be performed by children without supervision.

Children should be supervised to ensure that they do not play with the appliance.

## **SAVE THESE INSTRUCTIONS**

### **EQUIPMENT:**


- \* Irrigator
- \* 7 nozzles
- \* Nozzles box
- \* Travel bag
- \* USB cable
- \* Adapter
- \* User manual (warranty card is included)

### **APPLIANCE DESCRIPTION**

Model	POI-H350W	Input	5 B (1 A)
Material	ABS, PC	Battery capacity	2500 mAh
Working Power	21 W	Water Pressure Range	40-140 psi
Charging time	4-5 hours	Waterproof level	IPX7
Reservoir capacity	300 ml	Water flow method	Pulsate water

## FITTINGS

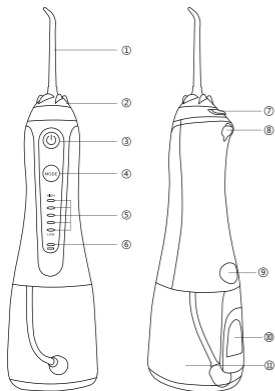- ① Nozzle tip
- ② Tip rotation dial
- ③ Power button
- ④ Mode selection button
- ⑤ Mode indicator lights
- ⑥ Charging indicator light
- ⑦ Tip release button
- ⑧ Charging port
- ⑨ Tank release button
- ⑩ Water tank cap
- ⑪ Water tank

## Various nozzles:

Tongue cleaner x1

Orthodontic x1

Periodontal x1

Classic x2


Toothbrush x1

For the nose x1


## CHARGING UNIT PRIOR TO FIRST USE

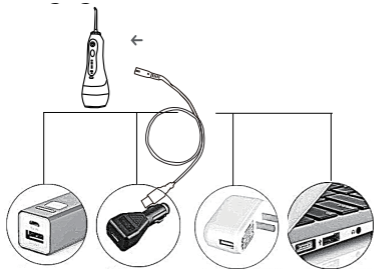This product comes with an AC/DC charger. Also the unit is charge through devices with 5V USB's power capabilities such as mobile phone, lap top, USB power wall charger, etc. The charge indicator on the handle will flash when charging and will stay on continuously when the unit is fully charged. A complete recharge will take 4-5 hours.


## FILLING THE RESERVOIR

Lift flip top on reservoir and fill the reservoir with lukewarm water.
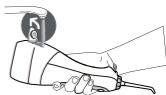
## INSERTING AND REMOVING TIPS

Insert tip into the center of the knob at the top of the Water Flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

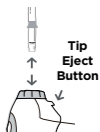**CAUTION: Do not press tip eject button while unit is running.**

## **ADJUSTING THE PRESSURE SETTING**

You can adjust the pressure setting on your Water Flosser by pressing the mode button on the handle. Start with the lowest setting and increase to your liking.

## **WHEN FINISHED**

Turn the unit off. Empty any liquid left in the reservoir.

## **TIP USAGE**

### **Jet Tip (Classic)**

Been over the sink and guide the jet to your teeth and gums, Glide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the sink. Direct the stream of water towards the gumline at 90 degrees, so that contact is made where the tooth and gums meet. Manoeuvre the flow of water from tooth to tooth remembering to pause briefly between each tooth.

### **Periodontal Tip**

Place the periodontal tip in contact between teeth and the gumline. Lean the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Start the irrigator by setting it at the lowest pressure and clean by gently following gumline, inserting the periodontal tip in between teeth. Continue tracing along the gumline to massage the gums and safely deliver solution into the gingival pockets below the gumline.

### **Orthodontic Tip**

Use the Orthodontic Tip to simultaneously brush

and rinse hard-to-reach areas around braces, crowns, bridges and implants. Start on the last molar (back tooth). Place tip at 90degree angle to the tooth at the gumline. Gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth and outside of both the upper and lower teeth. Continue process until you have cleaned the inside.

### **Tongue Cleaner Tip**

Lean your head over the sink. Glide the pressure button to the lowest setting. Start the irrigator and place the tongue cleaner at the back of your tongue and move it forward to clean your tongue. Repeat to cover entire tongue. Set the pressure to the lowest setting and adjust to a hinger setting for more water pressure as it feels comfortable.

### **Toothbrush Tip**

The hybrid device (irrigator with toothbrush tip) acts like an automatic toothbrush. It can be used with or without regular antibacterial paste. For the simultaneous brushing of teeth and irrigation. As part of complex measures at the first stage, a standard procedure for brushing teeth with a brush and antibacterial paste is provided. Only then is it possible to connect the irrigation system. At this stage, the alternate use of different types of attachments is allowed. The choice of models and the course of action is determined by the number and complexity of the dental problems present.

## **Standard Nasal Irrigation Tip**

Before the procedure for rinsing the nose, it is necessary to consult your ENT doctor for your specific case. There are contraindications.

Rinsing the nasopharynx with an irrigator takes less than 1 minute. It is convenient to combine it with daily brushing of your teeth. It is advisable to rinse the nose and immediately upon returning from places where there is a high risk of catching viral infections (bacterial environment).

To complete the procedure, you should:

- Set the pressure regulator of the irrigator jet to the minimum value - one and a half or two on the scale of division. Turn on the irrigator.
- Lean over the sink and turn your head to one side. To the nostril, which turned out to be higher, attach a nasal nozzle and turn on the solution supply on the irrigator handle. This will cause fluid to flow out of the lower nostril.
- Blow your nose after irrigating the first nostril. Then you need to change the position of the head and apply the solution to the other nostril.
- At the end, you need to blow your nose well to remove the remnants of the rinsing solution and mucus.
- Switch off the irrigator.

## TROUBLESHOOTING

Problem	Cause	Corrective action
Reservoir leaks	Seal ring is abnormal	Remove the seal ring and reinstall
Low water pressure	Low battery	Recharge the device
	Nozzle head blocks	Replace with new nozzle heads
Lasts only few minutes with full charge	Dead battery	Contact authorized service center
	Defective Charging cord	Try with new charging cable
	Adapter broken	Try with new adapter
Motor makes sound but no water coming out	Nozzle head blocks	Replace with new tips
	Suction valve dislocation	Keep device on, fill water to the hole where nozzle head goes into
Motor doesn't make sound	Low Battery (Keep the device idle for more than 3 month)	Activate the battery by recharging it for 16 hours
Product stopped working suddenly	Overheating protection is activated (continuous use 5 tanks of water with highest pressure)	Restart the device in 2 hours

**Note:** Excessive heat will cause machine to stop when in use, switch off machine, cool for 2 hours and re-start.

**If your problem still hasn't been fixed by above troubleshooting, please kindly contact with authorized service center for help.**

It may occur below situation when first time use this machine or left it unused for a long time: The nozzle head won't spray water immediately after turning on the machine. In this case, please follow below instructions to make it work:

Insert the tip in and press power button to turn on the device.

Add proper amount of water into the place where tips are inserted in.

## MAINTENANCE

Make sure to throw out the water inside the tank, and run the appliance for 2 to 3 seconds to drain the water after usage. Otherwise bacteria could collect in the remaining water. Do not submerge the unit in water or any other liquid. Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser. Reservoir may be washed with soapy water. Wash the dirt off from inside using a brush or similar, and wipe off with a cloth. Dry the reservoir when not using for a long period.

## BATTERY / PRODUCT CORRECT DISPOSAL

The unit contains a Lithium battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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## Warranty conditions:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual.

Service life of the product is 5 years. In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card.

If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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# WARRANTY CARD

# ARDESTO

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this particular instrument.

**We ask you to keep the coupon during the warranty period. When purchasing a product, require a full warranty card.**

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. \*
2. The service life of household appliances is 5 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
  - misuse and non-consumer use;
  - mechanical damage;
  - damage caused by the ingress of foreign objects, substances, liquids, insects;
  - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
  - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
  - when operating equipment in the power supply network with a missing single ground loop;
  - in case of violation of seals installed on the product;
  - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

\* Tear-off maintenance tickets are provided by an authorized service center.

**The completeness of the product is checked. I have read the terms of the warranty service, no complaints.**

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