# ARDESTO

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## IMPORTANT SAFEGUARDS

---

The Power toothbrush is safe and effective when used according to directions. Like any electrical appliance, however, basic safety precautions should be followed, especially when children are present.

### DANGERS-WATER

---

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not place or store the charger base where it can fall or be pulled into a bathtub, sink or toilet.
- Do not place or drop the charger into water or other liquid.
- Do not reach for a product that has fallen into water or other liquid. Unplug immediately.
- Do not use while bathing.
- Unplug charger base before cleaning and assure that it is dry before plugging it into an electrical outlet.

### DANGERS-VOLTAGE

---

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug the charger base into an outlet with a voltage other than specified on the bottom of the charger.
- This product is designed to operate within a range of 100-240 volts.
- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.

### DANGERS-WARNINGS

---

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This device can be used by children over 8 years of age and persons with physical or sensory or mental disabilities or with insufficient experience and knowledge, if they are under constant supervision or instructed in the safe use of the device and they understand the possible dangers.
- Children should be supervised to ensure they do not play with the product. Cleaning and customer service should not be performed by children without supervision.

- Use this product only for its intended use as described in this manual. This product is designed to clean your teeth, gums and tongue only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- The toothbrush is a personal care appliance and is not intended for use on multiple patients in a dental office or institution.
- If your toothpaste includes peroxide, baking soda or other bicarbonate (common in whitening toothpastes), thoroughly clean your brush head and handle to avoid cracking the plastic.
- Never operate this product if it has a damaged cord, plug or charger base.
- Do not use this product if it is not working properly, if it has been dropped or damaged or dropped into water. Immediately discontinue use if the product appears damaged in any way.
- Keep the toothbrush, charger base and power cord away from heated surfaces.
- If you experience discomfort or pain in your mouth, discontinue use of this product and contact a dentist.
- Do not attempt to remove the battery. The device has a built-in rechargeable battery that must be disposed of separately only after the end of the life of the electric toothbrush. More details in the section «Disposal».
- Do not clean the toothbrush in the dishwasher.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- Regarding how to clean the tooth brush surfaces in contact with hand and skin, please refer to section <MAINTENANCE>.
- The battery should be disposed safely.

### **Simplified Declaration of Conformity:**

Shenzhen Risun Technology Co., Ltd. hereby declares that the type of radio equipment (electric toothbrush) complies with the Technical Regulations for Radio Equipment. The full text of the declaration of conformity is available on the website at: <https://erc.ua/uk/declarations/>

frequency band(s) in which the radio equipment operates: 148.5 - 5000 kHz  
maximum radio-frequency power transmitted in the frequency band(s) in which the radio equipment operates: minus 15 dB $\mu$ A/m  
Software version: software does not apply

## MEDICAL-WARNINGS

---

- If you have had oral surgery or gum surgery within the past two months, consult your dental provider before using the product.
- Consult your dental provider if excessive bleeding occurs while using this product or minor bleeding continues to occur after one week of use.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).
- Consult your physician prior to using the toothbrush if you have any medical concerns.
- The toothbrush has been tested and is certified as per ETL and CE regulations.
- Always unplug this product immediately after using.
- Never use while sleeping or drowsy.

Production date and device serial number are indicated on the packaging.

## EQUIPMENT

---


Electric Toothbrush - 1 pc

Brush head - 2 pcs

Inductive charging base - 1 pc

Instruction manual (warranty card is included) - 1 pc

# CONFIGURATION


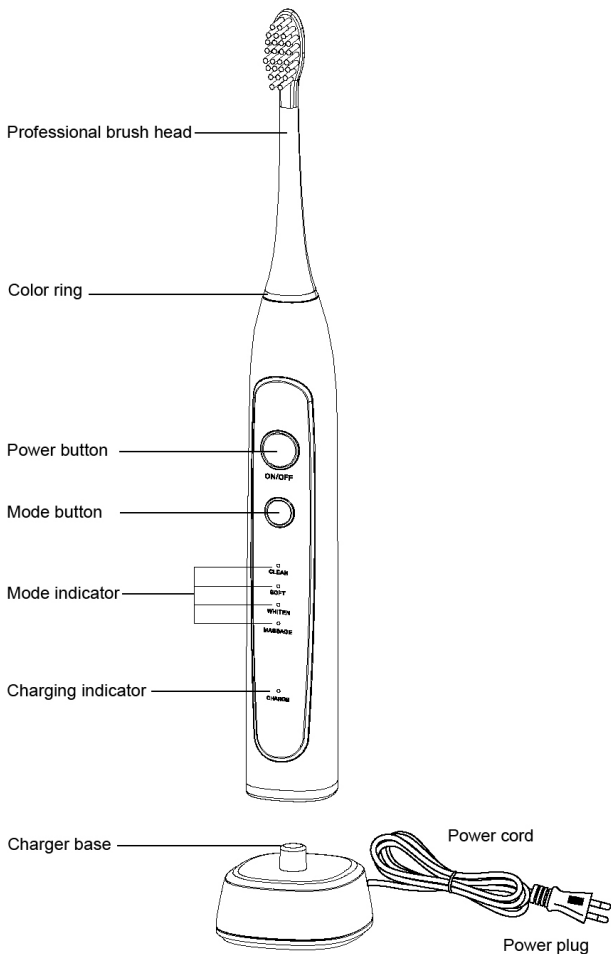
## GETTING STARTED

### Charging the toothbrush

Be sure to turn off the toothbrush before charging, plug the power plug into the outlet. Place the handle in the charger. The charger indicator light will be lit, indicating the handle is charging. After the battery is fully charged, charge indicator light will turn solid green, indicating the toothbrush is fully charged.

#### Notes:

- ① It takes at least 24 hours to charge the battery fully.
- ② The toothbrush must be charged for 24 hours, if:
  - a. this is the first time using the toothbrush;
  - b. the toothbrush has not been used for more than three months;
- ③ To keep the battery fully charged, we recommend that you keep the toothbrush in the charger when not in use.
- ④ During charging, the toothbrush will get warm, but this is normal.
- ⑤ The operating duration will become shorter as the battery gets older. If the toothbrush operates only for 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle. Contact service center for support.
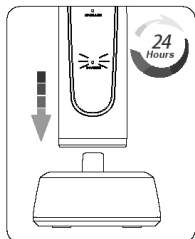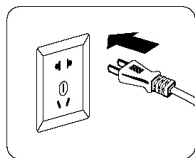
## HOW TO USE

- ① Slide the brush head onto the handle. Push brush head down until it locks into place on handle. Note: There is small gap between color ring and the handle. (Fig 1)
- ② Wet the bristles and apply a small amount of toothpaste.
- ③ Place the bristles against the teeth at a slight angle towards the gum line. (Fig 2)
- ④ Press the mode button to select mode.

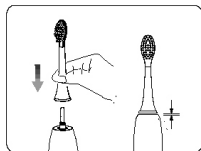
Fig. 1

If you know the mode used last time, you need not select mode again. There are 4 modes.

- a. **CLEAN:** Standard mode for superior teeth cleaning.
- b. **SOFT:** Gently caress, carefully clean sensitive mouth.
- c. **WHITEN:** Teeth cleaning with polishing and whitening.
- d. **MASSAGE:** Pulsates for gentle gum stimulation.

### Notes:

During brushing, if press Mode button, toothbrush will work at next mode and start another 2 minutes cycle and the mode indicator will light up.

- ⑤ Apply light pressure and gently move the brush head slowly across the teeth in a small back and forth motion.
- ⑥ In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section.
- ⑦ The toothbrush will automatically stop after 2 minutes of use. (Fig 3)
- ⑧ Pull brush head out from the toothbrush.(Fig 4)
- ⑨ Rinse both the brush head and handle under running water wipe them dry. (Fig 5)  
Put the handle back into the charger.

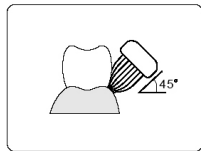
Fig. 2

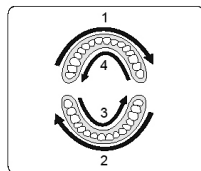
Fig. 3

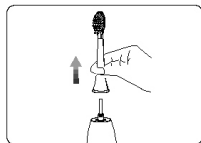
Fig. 4

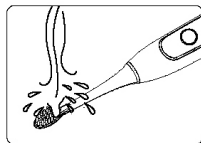
Fig. 5

## MAINTENANCE

Be sure to remove the plug from the electrical outlet before cleaning the product.

### Toothbrush

- ① Remove the brush head.
- ② Wipe all surface areas with a damp cloth and place the handle back into the charger.
- ③ Do not submerge the handle into water. Rinse the brush head and bristles after each use.
- ④ Do not clean in the dishwasher.

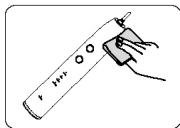
Fig. 6

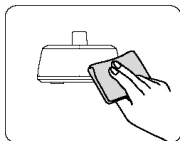
Fig. 7

### Charger

- ① Using a soft cloth, wipe off any dirt on the charger and the storage.
- ② Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause an electric shock, overheat or fire.


## TROUBLESHOOTING

Problem	Possible cause	Solution
The toothbrush doesn't work.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 24 hours.
Weak or diminished vibration.	The toothbrush is pressed hard against your teeth.	Place the brush head gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The battery has reached the end of its life cycle.	Contact service center for support.
	Charging time was insufficient.	Charge the main unit for at least 24 hours.


## DISPOSAL

The power source of this product is Lithium-ion battery which can be recycled. At the end of the product's life, always remove the battery before you discard the product, and recycle or dispose of the battery according to local waste management requirements. If you have recycling questions, contact your local waste management office.


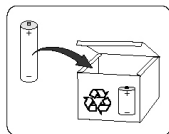
### Removing and recycling the battery

Do not disassemble the product at any time other than when disposing of it. Doing so earlier will destroy the product. Observe basic safety precautions. The appliance must be disconnected from the supply mains when removing the battery.

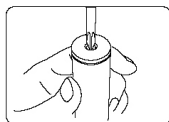1. Before removing the battery, be sure to completely drain the unit of power.
2. Use a screwdriver to remove the screw cover and battery door.
3. Hold the handle upside down, press the toothbrush shaft, and remove the two hooks to release the internal components of the handle.
4. Insert the screwdriver next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.
5. Recycle or dispose of the battery properly. Ask your local waste management office if they accept the remaining toothbrush components for recycling.

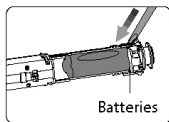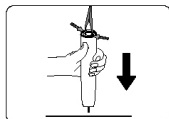Be careful not to short-circuit the negative and positive poles of the battery. Do not connect the appliance to the mains again after you have removed the battery.

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## Warranty conditions:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual.

Service life of the product is 3 years. In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card. If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by USC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card.

In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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# WARRANTY CARD

Dear Buyer! Congratulations on your purchase of the ARDESTO brand appliance, which was designed and manufactured in accordance with the highest quality standards, and we thank you for choosing this particular instrument.

We ask you to keep the coupon during the warranty period. When purchasing a product, require a full warranty card.

1. Warranty service is carried out only if there is a correctly and clearly filled original warranty card, which indicates: product model, date of sale, serial number, warranty service period, and the seller's seal. \*
2. The service life of household appliances is 5 years.
3. The product is intended for consumer use. When using the product in commercial activities, the seller / manufacturer does not bear warranty obligations, after-sales service is performed on a paid basis.
4. Warranty repair is carried out within the period specified in the warranty card for the product in an authorized service center on the conditions and terms determined by applicable law.
5. The product is withdrawn from the warranty in case of violation by the consumer of the operating rules set forth in the instruction manual.
6. The product is removed from warranty service in the following cases:
  - misuse and non-consumer use;
  - mechanical damage;
  - damage caused by the ingress of foreign objects, substances, liquids, insects;
  - damage caused by natural disasters (rain, wind, lightning, etc.), fire, domestic factors (excessive humidity, dust, aggressive environment, etc.)
  - damage caused by non-compliance of power and cable network parameters with state standards and other similar factors;
  - when operating equipment in the power supply network with a missing single ground loop;
  - in case of violation of seals installed on the product;
  - lack of serial number of the device, or inability to identify it.
7. The warranty does not cover consumables and accessories, as well as filters, shelves, drawers, stands and other containers for storing products.
8. The warranty period is 12 months from the date of sale.

\* Tear-off maintenance tickets are provided by an authorized service center.

The completeness of the product is checked. I have read the terms of the warranty service, no complaints.

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